


곰팡이 독소 8종

Cadenza CD-C18 HT 2.0 x 150 mm, 3 μm

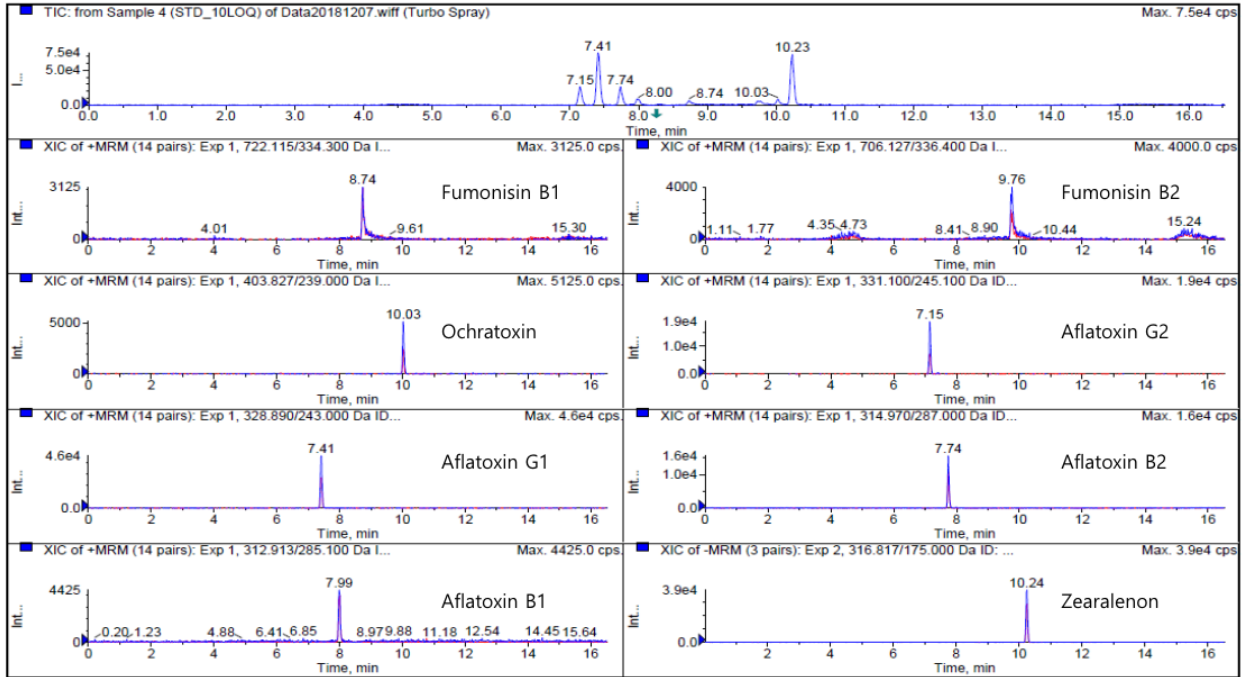

이동상 : A – 5 mM Ammonium formate / Formic acid = 100 / 0.1

B – 5 mM Ammonium formate in Methanol / Formic acid = 100 / 0.1, Gradient

온도 : 40 °C

검출기 : LC/MS-MS

주입량 : 5 μL

Time	Flow rate (μL/min)	B(%)	A(%)
Int.	200	5	95
0.5	200	5	95
2.0	200	40	60
9.0	200	100	0
12	200	100	0
12.1	200	5	95
16	200	5	95