


Desmosine, Isodesmosine

Intrada amino acid 3.0 x 100 mm, 3 μm

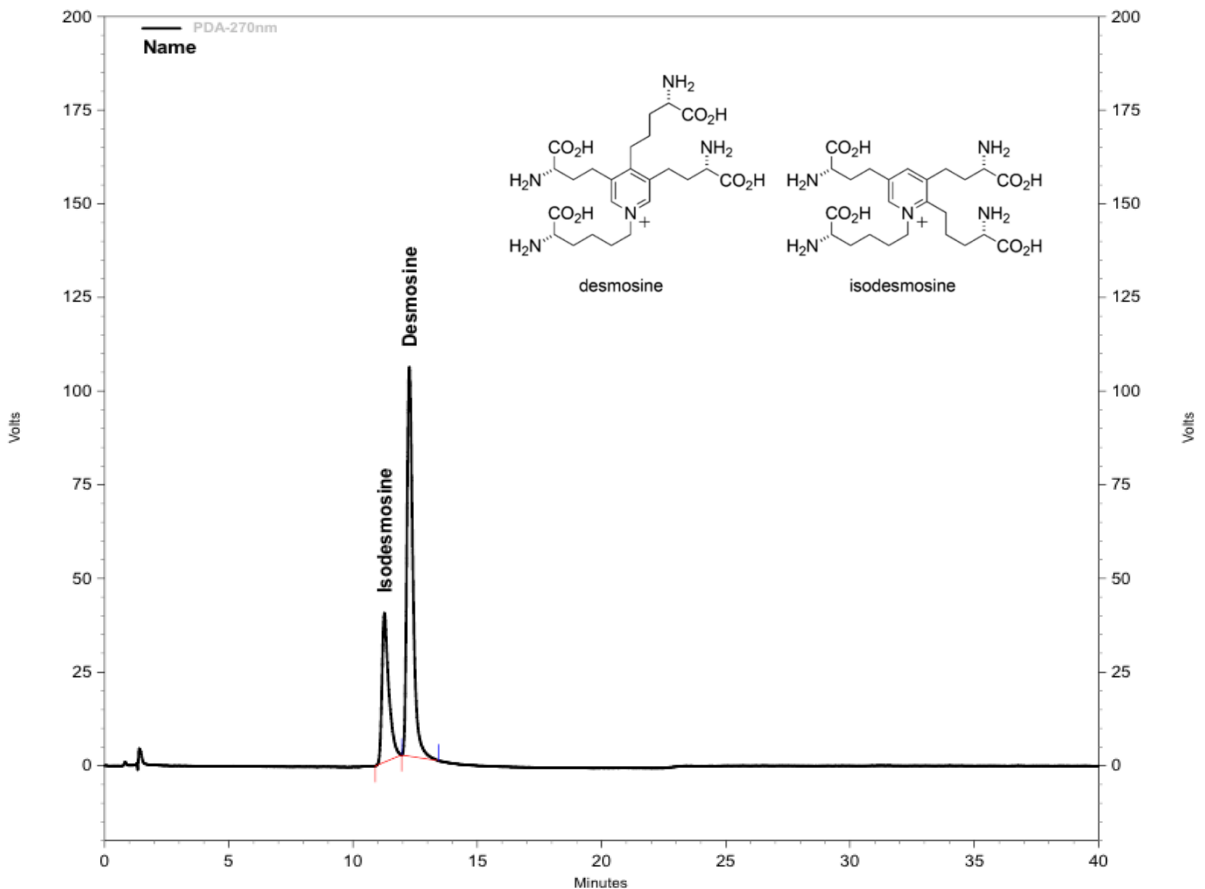

이동상 : A – 10 mM Ammonium formate

B – 100 mM Ammonium formate / Methanol = 70 / 30

온도 : 35 °C

검출기 : PDA (270 nm)

주입량 : 10 μL

샘플 : 150 μg/mL

Time (min)	Flow rate (μL/min)	B (%)	A (%)
Int.	400	20	80
15	400	50	50
20	400	50	50
21	400	20	80
40	400	20	80