


4 types of preservatives

보존료 4종

1. Isopropyl paraben (50 µg/mL)	2. Propyl paraben (50 µg/mL)
	
3. Isobutyl paraben (50 µg/mL)	4. Butyl paraben (50 µg/mL)
	

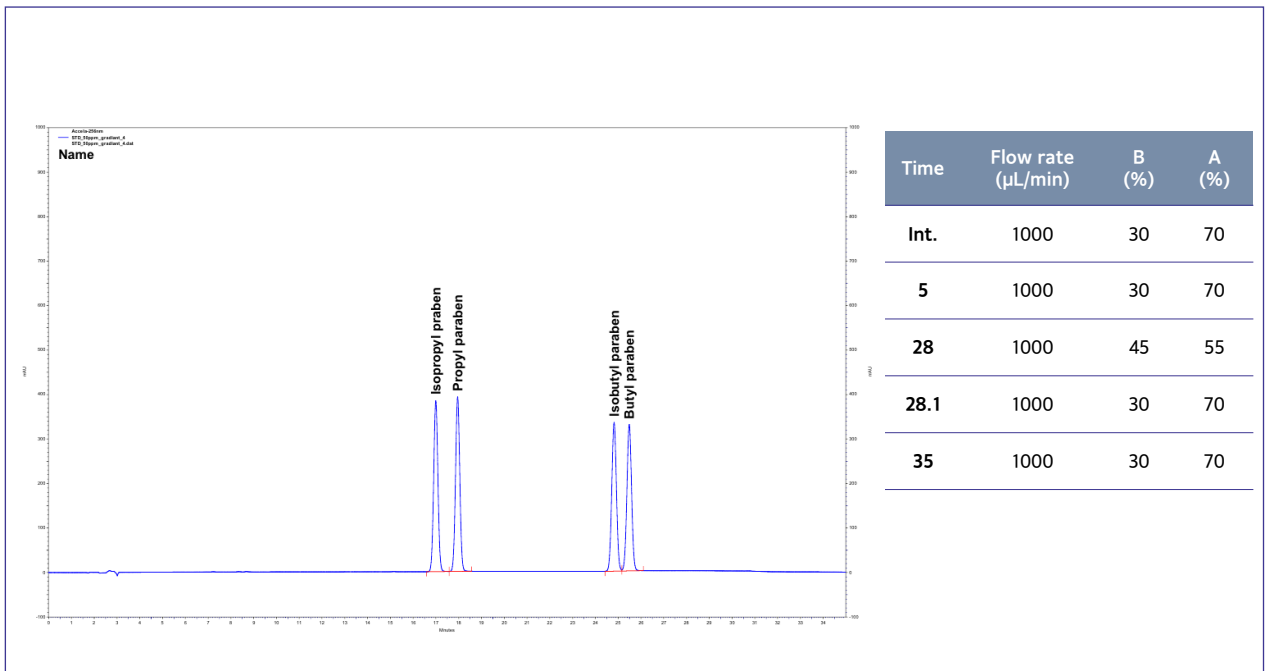

[HPLC Conditions]

System : NANOSPACE SI-2, OSAKA SODA, JAPAN
 Column : CAPCELL PAK C18 UG120 (4.6 ×250 mm, 5 µm), OSAKA SODA, JAPAN
 Mobile phase : A : DW
 B : ACN
 Flow rate : 1 mL/min
 Temperature : 35°C
 Detection : PDA (256 nm)
 Inj. vol. : 5 µL
 Standard : 제공받은 표준품을 Methanol 로 50 ppm 이 되게 희석하여 사용한다.