


4-Aminophenol

4-아미노페놀

Unison UK-C18

1. 4-Aminophenol
(10 µg/mL)

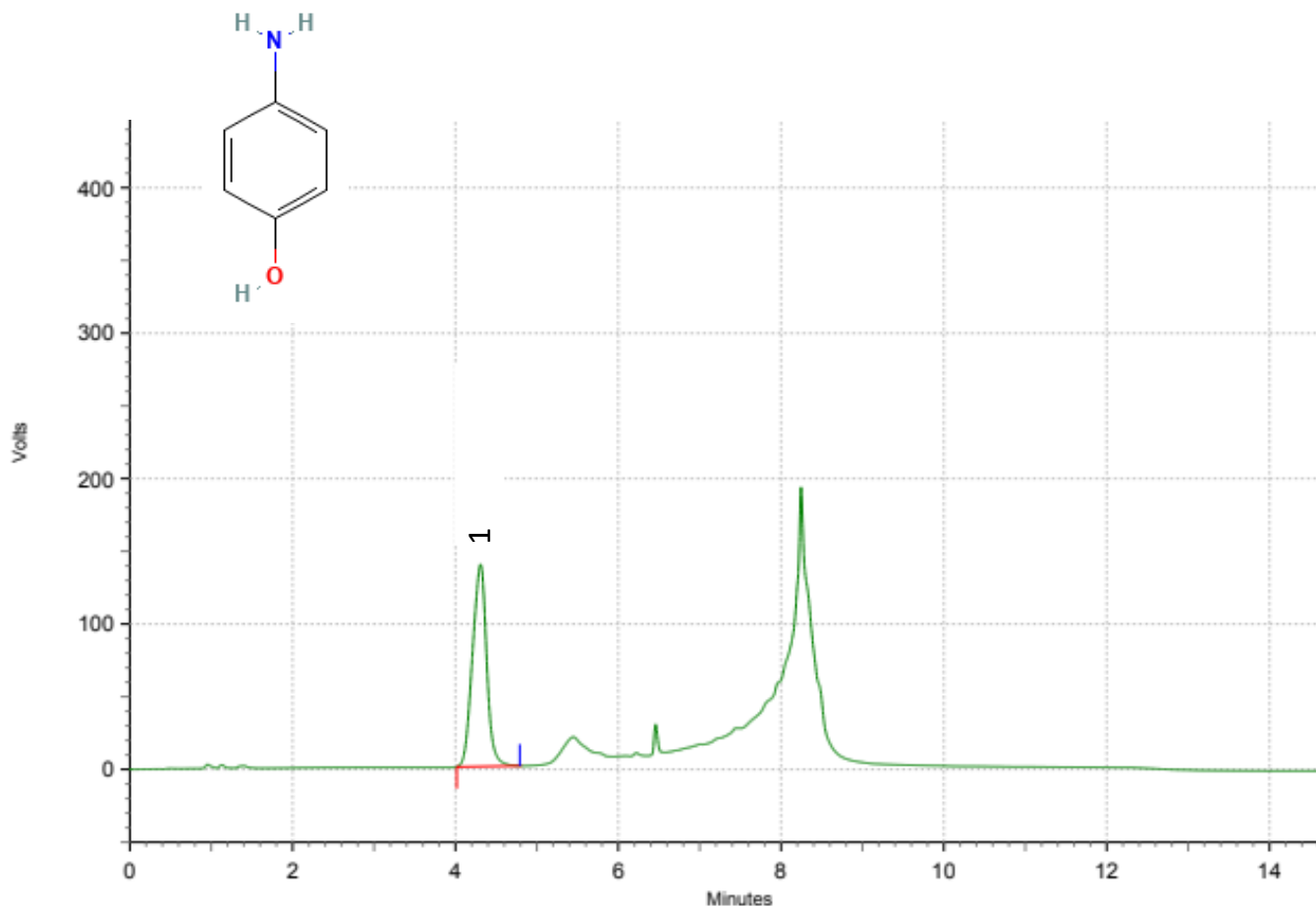

[HPLC Conditions]

Column	: Unison UK-C18 (4.6 × 100 mm, 3 µm), Imtakt, JAPAN
Mobile phase	: A : KH ₂ PO ₄ 1.7 g + Na ₂ HPO ₄ 1.8 g in DW B : MeOH 1%B (0-3min), 1-81%B (3-7min), 81-1%B (7-7.1min), 1%B (7.1-10min)
Flow rate	: 1.0 mL/min
Temperature	: 35°C
Detection	: PDA (230 nm)
Inj. vol.	: 10 µL (dissolved in MeOH)

