


Mirabegron

미라베그론

Cadenza CD-C18

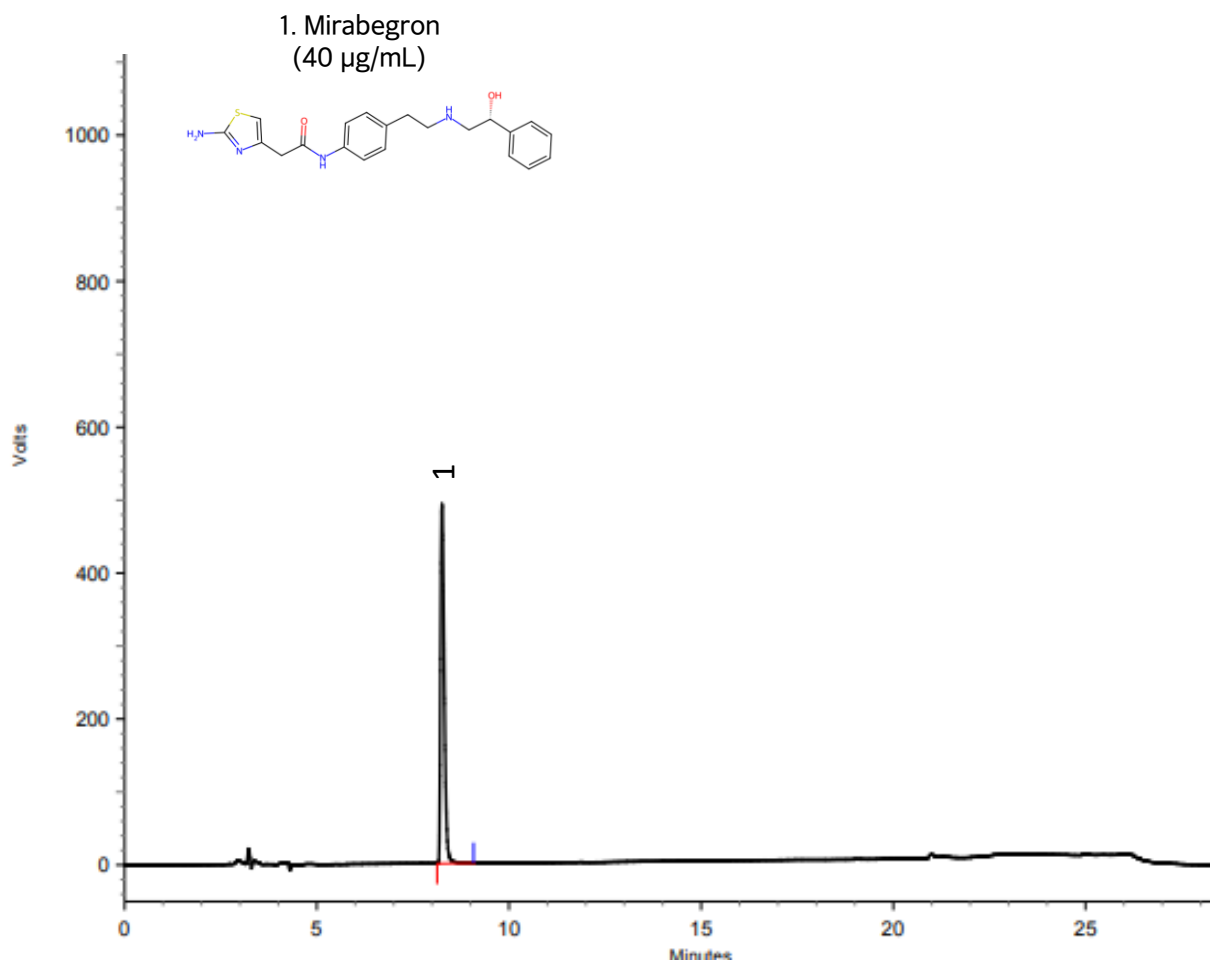

[HPLC Conditions]

Column	: Cadenza CD-C18 (4.6 × 250 mm, 3 µm), Imtakt, JAPAN
Mobile phase	: A : 10 mM KH ₂ PO ₄ (aq) B : ACN 20-70%B (0-23min), 70-20%B (23-23.1min), 20%B (23.1-30min)
Flow rate	: 0.8 mL/min
Temperature	: 35°C
Detection	: PDA (225 nm)
Inj. vol.	: 10 µL

